

Blood Collection Monitor BCM 20 Ultra




BCM 20 Ultra


REMI BCM-20 Ultra, Next generation blood collection monitor with interactive multicolor LCD display. Helps technicians to follow the SOP's ensuring standardization of the collection process

Removable magnetic tray ensures easy maintenance & transportation. The long lasting battery backup of 24 Hrs makes it an ideal choice for outdoor camps.

Self Diagnosis & System Checks

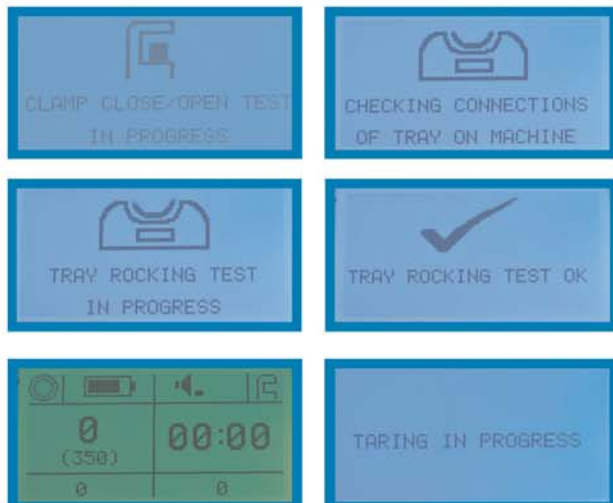


Interactive Multicolor Process Guide

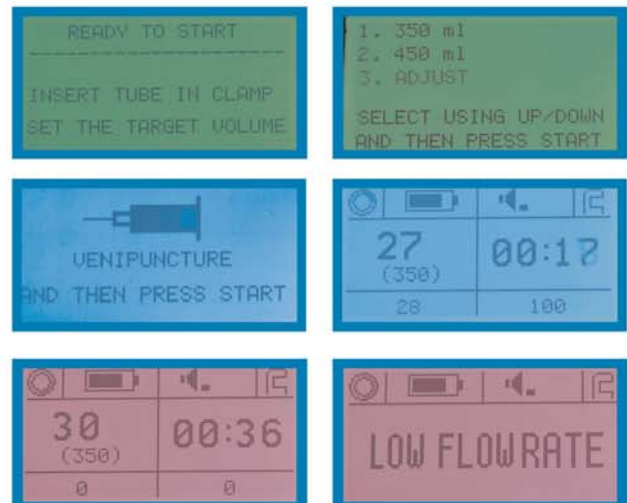

Blood Collection Monitor BCM 20 Ultra

Salient features

- Interactive multicolor display
- Adjustable alarm volume
- Encapsulated design allowing easy cleaning. No liquid intrusion in the device
- Carrying Case for device including all accessories
- Integrated Calibration function with the use of weights
- Pause facility to temporarily suspend the donation process
- Precise control on volume collection (Microprocessor based system), best in class accuracy (deviation $\leq 1\%$)
- Removable Magnetic Cradle: single grip tray removal
- Automatic release of the bag when lifted
- Audio Visual Alarm
 - Flow rate (< 20 ml/min or > 180 ml/min),
 - Completion of collection
 - Battery depletion
- Visual Indication & Displays
 - Set volume (ml)
 - Actual volume (ml)
 - Collection Process status (Run / Stop / Pause)
 - Actual weight (gm)
 - Flow Rate (ml/min)
 - Time of collection (min)
 - Battery status

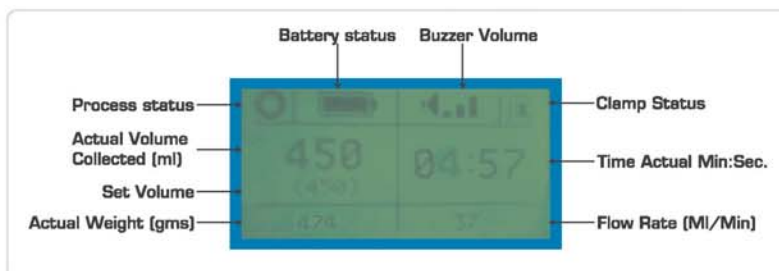

Technical Specifications

Parameter		BCM 20 Ultra
Mixing	Cycle / minute	12-14
Battery back up/ rating	Hours	Up to 24 (with overcharge protection)
Max. Volume	ml	Up to 600 (net collection)
Taring range	gm	Up to 600
W x D x H	mm	275 x 318 x 205
Net Weight	Kg.	5
Electrical	220-240 volts, 50 Hz, SINGLE PHASE	

Specifications subject to change without notice

V. S. International Pvt. Ltd.

A-204, Neelam Centre, Hind Cycle Road, Worli, Mumbai - 400 030., INDIA.
Contact : +91 22 4333 3900 / 2498 8777 Fax: +91 22 6660 8779